



Enzymax LLC, 870 Corporate Dr., Suite 201; Web: www.enzymax.net; Tel: 859-219-8482; Fax: 859-219-0653

Description:


EZC222F Midi Filter columns are designed to work with DNA or RNA midi columns. It is used for clarification of lysate during DNA/RNA midi preparation with midi DNA/RNA binding columns (either silica based or ion-exchange column).

Specifications:

- Polypropylene housing
- Used with vacuum manifold
- Column capacity: 14 ml
- Binding: do not bind to DNA, RNA
- Filter membrane pore size: 90 μm
- Can be attached to DNA midi binding column (see picture below) to filter the cell debris and clarify the lysate, which prevents the clogging of DNA midi binding columns.

Pictures:

Price: \$30 for pack of 25 columns (Cat#EZC222F)

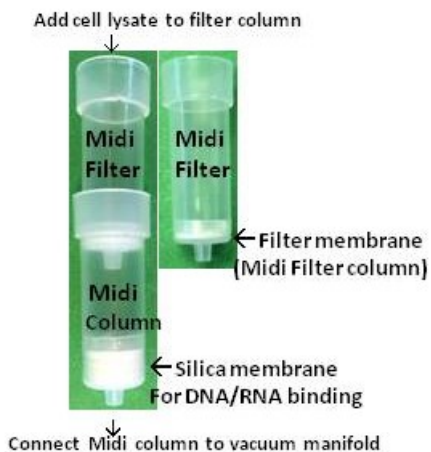

Related Product: DNA and RNA Midi spin column (EZC222and EZCR222)

http://www.enzymax.net/columns/midi_column_DNA_RNA.htm

DNA/RNA Midi Spin Column

